


# CFPMFXH

IN STOCK PRODUCT

Visit [www.oakworksmed.com](http://www.oakworksmed.com) for pricing.

- **Best value in full carbon fiber top design**
- **500 lb. Patient Weight Capacity**
- **56" of Metal Free Imaging Area**
- **4 Year Warranty**

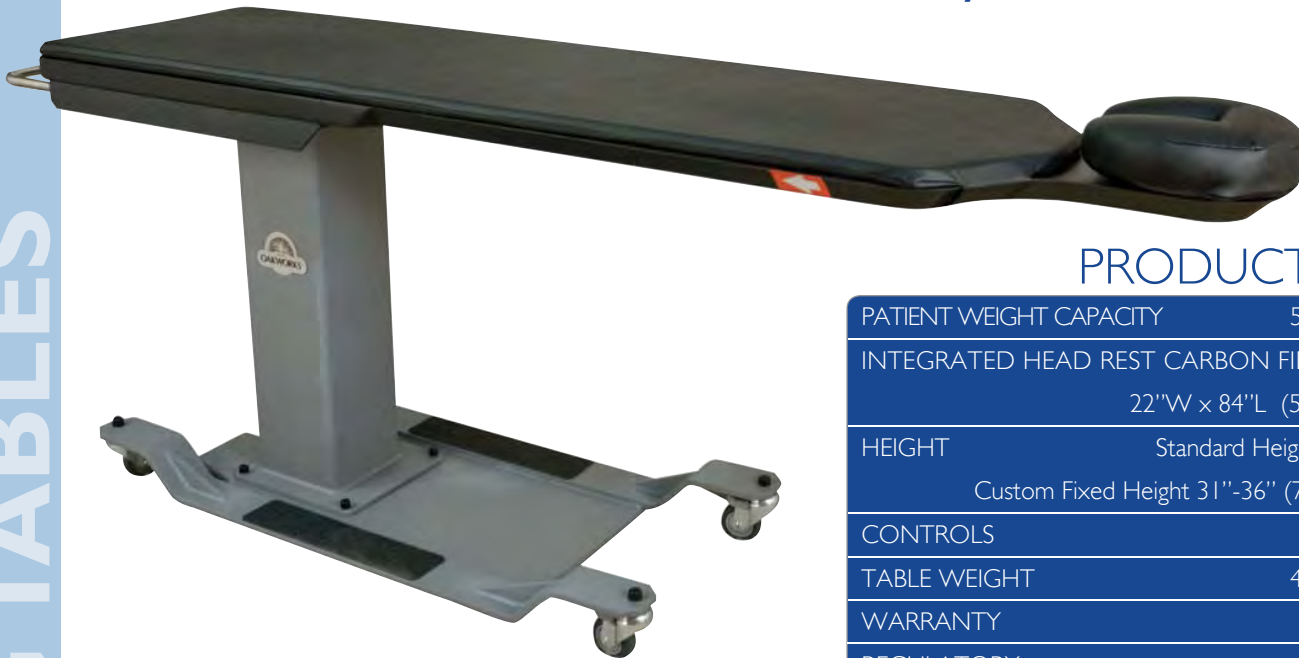

## PRODUCT SPECS

PATIENT WEIGHT CAPACITY	500lbs. (249kg.)
INTEGRATED HEAD REST CARBON FIBER TOP	22"W x 84"L (56cm x 213cm)
HEIGHT	Standard Height 37" (94 cm) Custom Fixed Height 31"-36" (79 cm - 91 cm)
CONTROLS	n/a
TABLE WEIGHT	415lbs. (188kg.)
WARRANTY	4 year
REGULATORY	FDA, ETL, CE

The New Oakworks® CFPMFXH Fixed Height C-Arm Imaging Table incorporates Oakworks brand reputation for high quality and reliability into this cost effective table designed to accommodate mobile and ceiling suspended C-Arm systems. The CFPMFXH is offered in a standard height of 37" with custom height options available in one inch increments from 31" to 36". This mobile c-arm imaging table has a 500lbs patient weight capacity and provides all of the basic features required for successful pain management image-guided procedures. The integrated head rest carbon fiber top features a 56"(142cm) metal free imaging area for maximum C-Arm efficiency. The table is mounted on a low profile base and can be moved easily from room to room and secured with 4 locking casters. Included with every table are a 1" radiolucent table top pad (2" pad is optional), 1 crescent face rest pad, and 2 supine cranial pads. See page 1 for pad descriptions.



HAND RAIL  
Move your table with ease

## RECTANGULAR TOP OPTION

The rectangular top is wider and shorter (24"x80", 61cm-203 cm) than the integrated headrest top providing greater flexibility.

*The three cranial positioning pads are not included with this option.*



**OAKWORKS®**  
**MEDICAL**

Made in the USA   
with US & imported parts